



PARAMOUNT BED CO., LTD.
2-14-5 Higashisuna Koto-Ku, Tokyo 136-8670, Japan
<http://www.paramount.co.jp/english/>

PARAMOUNT BED INDIA PVT. LTD.
Plot #430, Sector-8, IMT Manesar-122050, Distt. Gurugram (Haryana)
Tel: +91-124-4205748



A5 series
Manual Bed

for a safer
healing
environment



PARAMOUNT BED INDIA commits to produce and deliver state-of-the-art, safe and reliable products designed in Japan through integrating Quality Management Principles and our experience of many decades.

Designs and specifications are subject to change without notice. All images of the product, its components and accessories may not reflect the final product. The images may reflect optional equipment and are subject to change without notice. All information presented herein is based on current data or design goals available at the time of posting and is subject to change without notice. Some equipment and/or features may be unavailable. Actual results of performance and other specifications may differ or vary with production models and may depend on selected models and/or accessories.

Product Features Highlight

Safety in Design

- The Paramount Bed A5 Series complies with the requirements stipulated in IEC60601-2-52 meant to reduce the risk of patient entrapment in and around side-rails.

Split Type Side-rails

- A low bed height of 350mm*1 ensures that patients can firmly place their heels on the floor while the ergonomic side-rail design provides patients with additional support when standing up. The beds feature a gapless design between the bed frame and side-rails for maximum stability as patients sit by the bed and when they stand up.

*1 Height from floor to mattress base surface



Anti-tampering Lock

- This safety feature keeps the patient safely within the bed by preventing accidental release of the side-rails while the patient is leaning against them. This feature also prevents the patient from falling over when he is attempting to lower the side-rail himself.

Angle Indicator

- Medical personnel can gather clear and easy readings of the back raise angle of the patient. The angle indicator assists medical personnel by ensuring that the backs of patients have been raised to above 30°, which prevents ventilator-associated pneumonia (VAP).




Easy Cleaning & Maintenance

- Maintenance staff and medical personnel can easily remove or raise sections of the mattress base for easy cleaning and convenient top-down parts replacement.



SSR (Stress-Saver Raise)™

- Stress on the patient's body is reduced during back raise, as the back section of the bed frame extends to accommodate a natural sitting up process. This also prevents the patient from sliding down towards the foot of the bed.

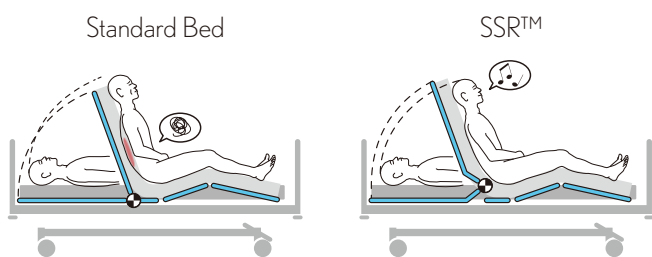

Collapsible Type Side-rails

- Visibility through this side-rail provides patients with a sense of openness while allowing medical personnel to be more readily aware of a patient's condition.




Side-rail position when lowered.



Side-rail position when raised.

Compact Design

- Other than providing a comfortable space for patients to rest on, the compact design of the overall bed width facilitates greater ease when manoeuvring the bed around sharp corners, narrow entrances and elevators. As part of its compact design, gaps between the bed and its components are reduced. This allows medical personnel to closely position a stretcher next to the bed for quick and easy transfer of patients.

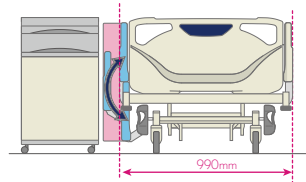

Casters



Caster with Individual Locking System Type



Leg Cap Type

[Individual Locking System]

- Depressing the stoppers located on each of the casters at the four corners allow you to lock down the bed.

Model No.		IA-57380	IA-58380	IA-59380
Type		1Crank	2 Cranks	3 Cranks
Total Length		2,200mm		
Total Width		990mm		
Total Height		955mm (with caster) / 885mm(with leg cap)		755 to 1,105mm
Mattress Base Height*1		550mm (with caster) /480mm(with leg cap)		350 to 700mm
Mattress Size		2,000 x 860 (mm)		
Safe Working Load*2		180kgf		
Back Raise Angle		0° to 70°		
Knee Raise Angle		-	0° to 25°	
Operation Control		Crank Handle		
Side-rail	Split	-	-	△
	Collapsible	○	○	△
Caster	Plastic	-	-	○
	Metal	△	△	-
Leg Cap		△	△	-

○ in the table is for standard feature. △ is for optional.

*1 Indicates the height from floor to the mattress base surface.

*2 Maximum load with which the bed can be operated.

(The total load of user's weight and incidental items including the mattress and options.)



Optional Items

Model No.	IZ-E001
Total Length	2,000mm
Total Width	860mm
Total Thickness	100mm
Materials	Cover: Polyester/ Padding: Urethane foam



PF-1200 Series



PF-3100A



KF-282