

SKANREVIVE

Biphasic Defibrillator



Salient Features

- Current controlled biphasic defibrillation.
- Synchronised cardioversion
- More than 100 shocks from a new fully charged battery.
- Charge to 300J in less than 8 seconds from mains or battery.
- Quick energy selection via rotary switch.
- Suitable for adult & paediatric applications.
- 6.5" colour LED display
- 3L/5L ECG support

CE 0843



RADIOLOGY



MONITORING



CARDIOLOGY



ANAESTHESIA



VENTILATION



DENTAL



SURGERY



TELE SOLUTIONS



VETERINARY

SKANREVIVE

Biphasic Defibrillator

Defibrillation technique

SKAN REVIVE employs current controlled, time limited biphasic defibrillation with patient impedance compensation to give accurate results.

- Reduces tissue damage by accurately controlling the current.
- Accurate delivery of energy.
- Biphasic pulse waveform with better efficiency compared with monophasic pulses.

Capable of delivering up to 300J, SKAN REVIVE is a powerful and reliable device available for reviving victims of SCA.



Efficiency

Thanks to the intelligently designed power module, SKAN REVIVE can charge to 300J in under 8 seconds when running from battery or AC mains. With a high capacity internal power source, SKAN REVIVE can deliver more than 100 full energy (300J) shocks from a new fully charged battery.

Quick and Easy

The ergonomically designed energy selector and the simple to use 1 2 3 operation helps to keep the focus where it is needed. With the rotary energy selector and illuminated energy indicator ring, setting the right energy gets quick, easy and accurate.




Patient monitoring

SKAN REVIVE is equipped with a 3 lead / 5 lead ECG for monitoring the patient after defibrillation therapy. SKAN REVIVE also has optional NiBP and SpO2 monitoring capabilities. With a host of monitoring features, SKAN REVIVE simplifies the task of post shock monitoring.

Intuitive User Interface

The large clutter free colour LED display on SKAN REVIVE helps users to quickly identify information provided on screen. The ergonomically placed keys with back-illumination makes operating SKAN REVIVE still simpler.

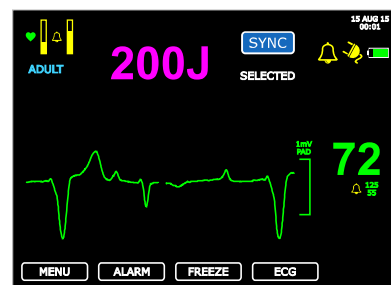

CUSTOMER INTERACTION CENTER

Toll Free No.: **1800-425-7002**
Email : cic@skanray.com

Monday to Saturday - 10AM to 6:45PM

Distributor:

Skanray Technologies Pvt Ltd.

Regd Office & Technology Division: #15-17 Hebbal Indl Area, Mysore - 570016, Karnataka, INDIA. Ph: +91 821 2415559, Fax: +91 821 2403344.
Healthcare Division: Plot #360, KIADB Industrial Area, Hebbal, Mysore - 570018, Karnataka, India. Ph: +91 821 2407000, Fax: +91 821 2407001.
Skanray Europe S.R.L: Via Cicogna 36 & 38, 40068 San Lazzaro di Savena, Bologna, ITALY. Phone: +39 051 6255019.
Email: enquiry@skanray.com | www.skanray.com | CIN: U72200KA2007PTC041774

Warning: The equipment is to be used by only qualified Medical Personnel and in accordance with instructions of the manufacturer, as specified in the User Manual. It should be serviced and maintained by the Company's authorized service personnel. Skanray continuously upgrades its technology, performance and safety features. Pictures and specifications featured in this catalog are only indicative. Please refer to the commercial offer document for exact specifications and performance.

